


Meet your new mPERS

Belle+

- 1 Speaker**
Used to hear the Care Center
- 2 Help Button**
Press to send emergency alarms to the Care Center
- 3 Battery Light**
The light will blink red when device is in the charging cradle. Once completely charged, the red light will remain solid
- 4 Call Light**
The light will turn blue when making a call
- 5 Wifi Light**
Indicates WiFi connection
- 6 GPS Light**
Will blink when the device is communicating to the GPS satellite
- 7 Microphone**
Used to communicate with the Care Center

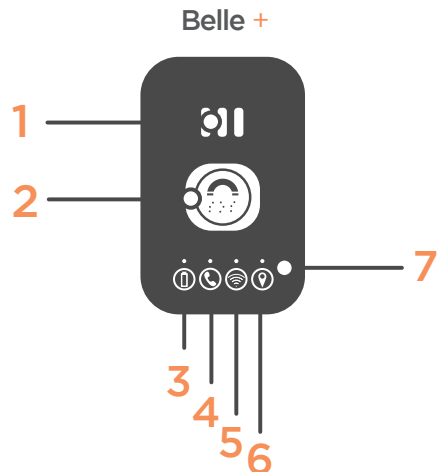

1400 Commerce Center Drive
Franklin, OH 45005

V03042019

If you know someone who could use our assistance, call us at **1-800-860-4230** or log on to **www.vricares.com** for more information on our products and services.

Belle+


Instructions

mobile medical alert systems




Belle + Instruction Guide

- 1 Plug the mPERS charging cradle into a power outlet. Make sure that the outlet is **NOT** operated by a light switch or a power strip. A red light on the cradle will light up.

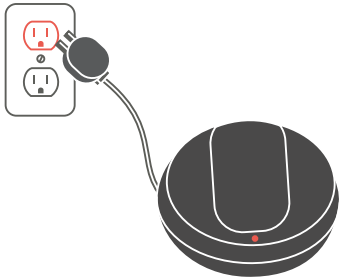

- 2 Place the mPERS pendant into the charging cradle. The red battery light will flash as the unit charges. The red light on the unit will be solid when fully charged.

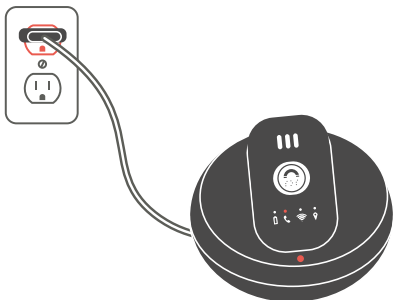

- 3 Once charged, put on your pendant. **Press and hold** your help button until the blue call light comes on. Your mPERS unit is now calling the Care Center.



- 4 The Care Center will talk to you through your unit. Let us know you are a new client testing your unit!



You're going places. Take us with you.

Belle+ is a lightweight, two-way communication pendant with a long lasting battery that gives you up to 3-5 days of coverage in one charge!



It is important to **test your device monthly** to ensure device is working properly. Before you test your unit, make sure it is **completely charged**.



Fall detection is not a perfect science. If you are able, **always** press your help button in the case of an emergency.

Belle+ Light Guide

Light	Icon	State	Description
Battery Light	((●))	Flashing	Flashing 1 second - low battery Flashing 2 seconds - charging
		Solid	Fully Charged
WiFi Light	((●))	Flashing	Device is acquiring WiFi signal when placing a call
		Off	Normal Operation
Call Light	((●))	Flashing	Normal Operation
		Failed	If light is solid for 5 seconds, then off for 5 seconds call has failed
GPS	((●))	Flashing	Device is acquiring signal when activated
		Off	After 5 minutes of being removed from cradle will turn off (normal)