


Quick Install

4G In Home Cellular Emergency Response System

Welcome to your new personal emergency response system (PERS).

If you have any questions or in need of assistance, do not hesitate to call us at [1-800-860-4230](tel:1-800-860-4230).


- 1 Plug the power cord into a power outlet. Make sure the outlet is **NOT** operated by a light switch or a power strip.

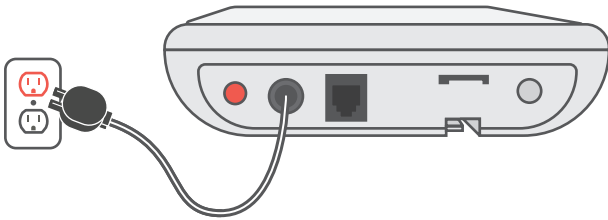

- 2 Press the red power button then wait for the unit to announce "System Ready".



- 3 Locate and put on your pendant or fall sensor, then press the button to activate your device. The unit will say "Calling The Care Center".

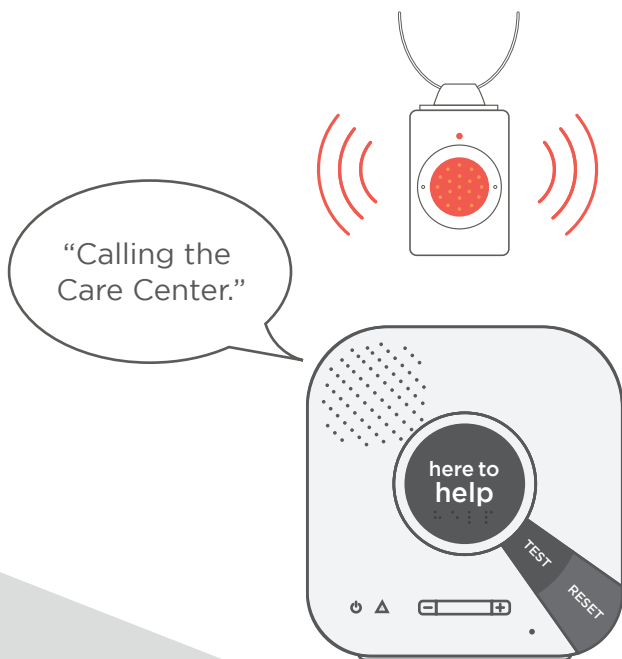

- 4 The Care Center will talk to you from your PERS unit. Let us know you are a new client testing your unit!



✦ If the Care Center does NOT respond within 2 minutes, call [1-800-860-4230](tel:1-800-860-4230) for technical assistance.