


Quick Install


CL4G Emergency Response System

Welcome to your new personal emergency response system (PERS).
If you have any questions or in need of assistance, do not hesitate to call us at [1-800-860-4230](tel:1-800-860-4230).

- 1 Plug the power cord into a power outlet. Make sure the outlet is **NOT** operated by a light switch or a power strip.

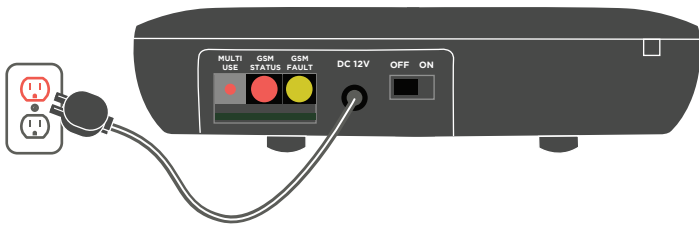

- 2 Switch the ON/OFF switch to the **ON** position and wait for the red GSM Status light to start blinking red every 2 seconds.

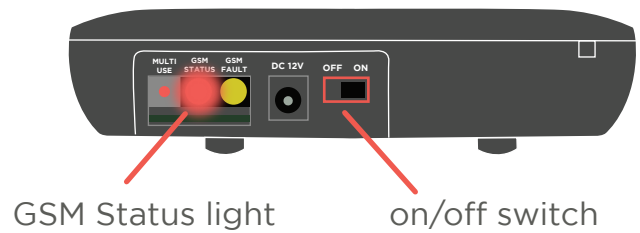

- 3 Locate and put on your CL4G pendant or fall sensor, then press the button to activate your device. The unit will say "Emergency call in progress".



- 4 The Care Center will talk to you from your PERS unit. Let us know you are a new client testing your unit!



✦ If the Care Center does NOT respond within 2 minutes, call [1-800-860-4230](tel:1-800-860-4230) for technical assistance.