


Meet Your New PERS

DCP100 Emergency Response System

DCP100 was created to provide you comfort, reassurance and peace of mind by making it possible for help to be just a touch of a button away. Simply connect your unit through your home telephone to begin living confidently knowing we are always here to help.

Get to Know Your Unit

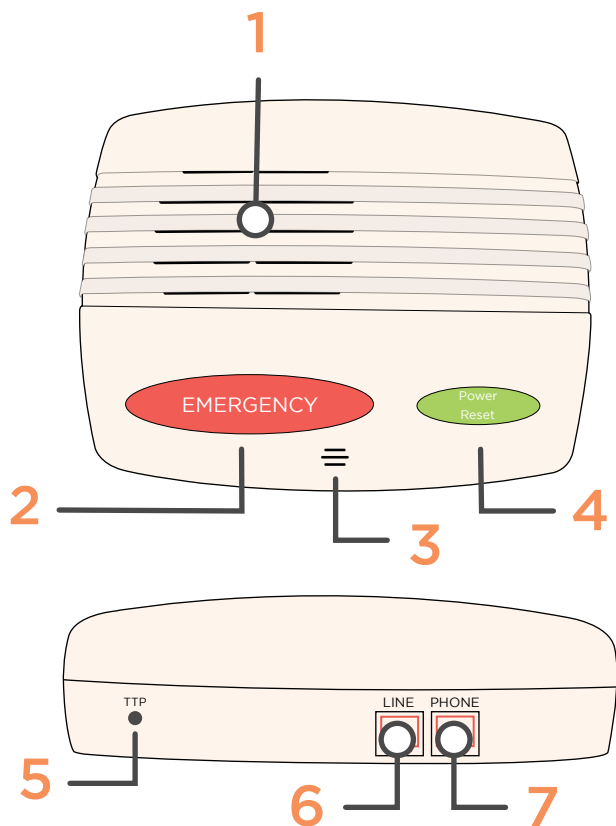

- 1 Speaker**
Used to hear the Care Center operators
- 2 Emergency Button**
Press to send emergency alarms to the Care Center
- 3 Microphone**
Used to talk to the Care Center operators
- 4 Power Reset**
Used to power off the DCP100 after power transformer has been unplugged from wall outlet
- 5 TTP**
For service providers use only
- 6 Line Jack**
Connects the unit to your telephone wall jack
- 7 Phone Jack**
Connects your unit to your telephone

Unit Facts & Features:

- ✦ The range on the PHB reaches up to 600 feet
- ✦ The button can be worn around your wrist or neck like a necklace
- ✦ Please test your unit monthly to ensure it is working properly
- ✦ Water-resistant button
- ✦ When connecting a "cordless phone" to the DCP100, make sure that the base for the cordless phone is a distance of 1 foot or greater between the base and the DCP100. Placing a cordless phone base too close to the DCP100 will cause a loud buzzing sound
- ✦ For additional questions, call tech support at [1-800-860-4230](tel:1-800-860-4230)