


Quick Install

DCP100 Emergency Response System


Welcome to your new personal emergency response system (PERS).

If you have any questions or in need of assistance, do not hesitate to call us at [1-800-860-4230](tel:1-800-860-4230).

- 1 Plug the power cord from the unit into an electrical outlet. Make sure the outlet is **NOT** operated by a light switch or power strip.

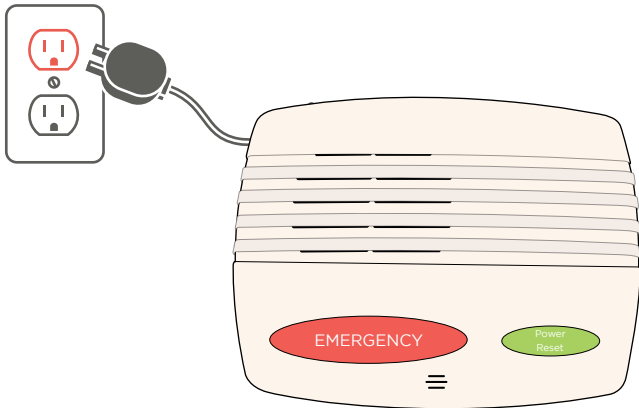

- 2 Connect the provided telephone cord into the DCP100 PERS jack labeled **WALL**. Insert the other end of the telephone cord into your home telephone wall jack.

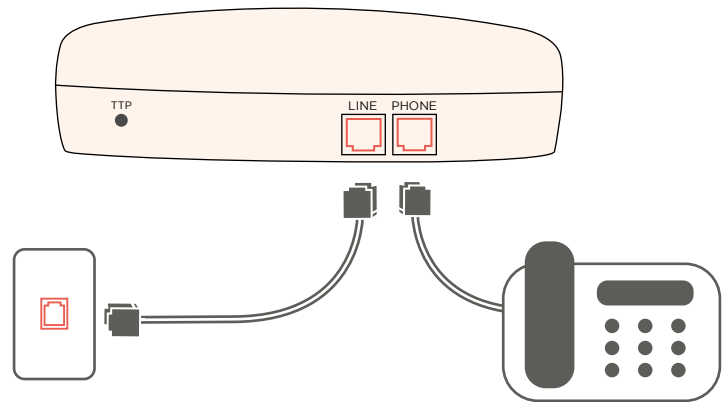

- 3 Locate and put on your DCP100 PERS pendant. **Press** your help button to activate your PERS unit.



- 4 The Care Center will talk to you from your PERS unit. Let us know you are a new client testing your unit!



If the Care Center does NOT respond within 2 minutes, call [1-800-860-4230](tel:1-800-860-4230) for technical assistance.