


Quick Install

Digital Emergency Response System

Welcome to your new personal emergency response system (PERS).

If you have any questions or in need of assistance, do not hesitate to call us at [1-800-860-4230](tel:1-800-860-4230).

- 1 Plug the power cord into a power outlet. If your home telephone is already plugged into a telephone outlet, remove the telephone cord from the outlet to insert in the unit **PHONE** jack. Then, connect the provided telephone cord into the unit **WALL** jack and your home telephone jack.

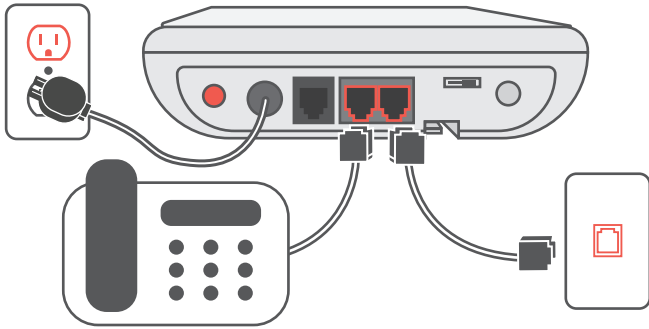

Make sure the outlet is NOT operated by a light switch or a power strip.

- 3 Locate and put on your pendant or fall sensor, then press the button to activate your device. The unit will say "Calling The Care Center".



- 2 Press the power button to turn on the unit. The device will announce "Power On".



- 4 The Care Center will talk to you from your PERS unit. Let us know you are a new client testing your unit!



✦ If the Care Center does NOT respond within 2 minutes, call [1-800-860-4230](tel:1-800-860-4230) for technical assistance.