


# Meet Your New PERS

## Linear 2400 Emergency Response System

Linear 2400 was created to provide you comfort, reassurance and peace of mind by making it possible for help to be just a touch of a button away. The Linear2400 was designed to be easy to install and maintain.

### Get to Know Your Unit

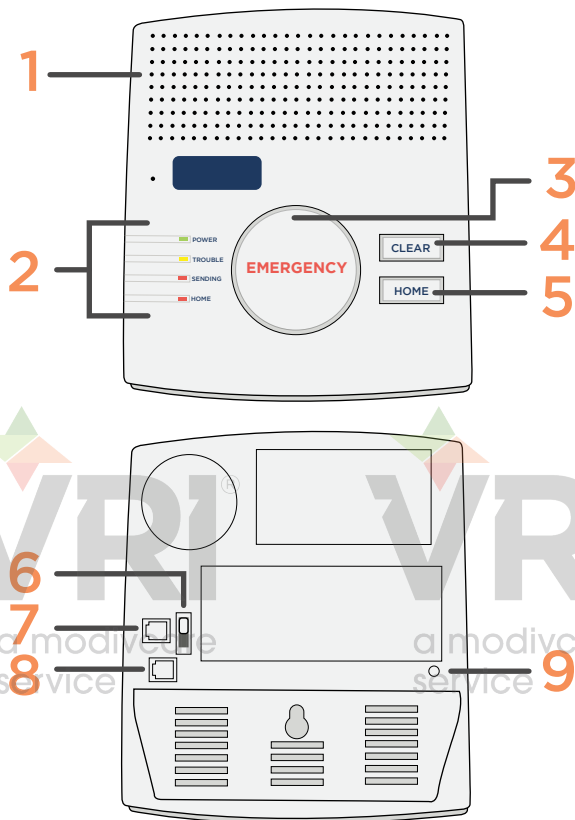

- 1 Speaker/Mic**  
Used to talk to the Care Center
- 2 Status Lights**  
Displays the current status of the unit
- 3 Help Button**  
Press to call the Care Center (works just like the button)
- 4 Clear Button**  
Press to stop alarms
- 5 Home Button**  
Press to reset activity timer
- 6 Power Switch**  
Flip to turn your unit on or off
- 7 Wall Jack**  
Used to connect your wall jack, or router
- 8 Phone Jack**  
Used to connect your home phone to the unit
- 9 Range Test Button**  
Used to test the button range around the residence

### Unit Facts & Features:

- ◆ The 2400 is compatible with standard phone lines
- ◆ The button can be worn around your wrist or neck like a necklace
- ◆ The button has a range of up to 600 feet, (average 3-5 year battery life), and is water resistant
- ◆ Care Center monitors unit's battery status
- ◆ 32+ hours of back up battery in case pf AC power loss
- ◆ **Status Lights by Color**
  - **Power Solid Green** - Unit has household power
  - **Power Flashing Green (Slowly)** - backup battery is charging
  - **Power Flashing Green (Quickly)** - AC power is not available, unit is running on backup battery
  - **Home (ON)** - Lights when activity timer is on
  - **Home Red (Blinks)** - When activity timer has less than 15 minutes left or has expired
  - **Sending Red (Flashing)** - Unit is trying to call the Care Center
  - **Trouble Yellow** - Low battery detected on your button

For additional questions, call tech support at 1-800-860-4230