
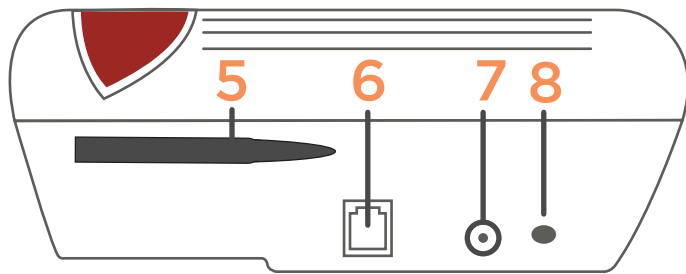


Meet Your New PERS

MXD-3G Emergency Response System

MXD-3G was created to provide you comfort, reassurance & peace of mind by making it possible for help to be just a touch of a button away. The MXD-3G was designed to be easy to install and easy to maintain. The unit has its own internal modem, so the MXD-3G does not use the clients cell plan in any way. No phone lines, no confusing connections, just plug it in household power and you are ready to go.

Get to Know Your Unit

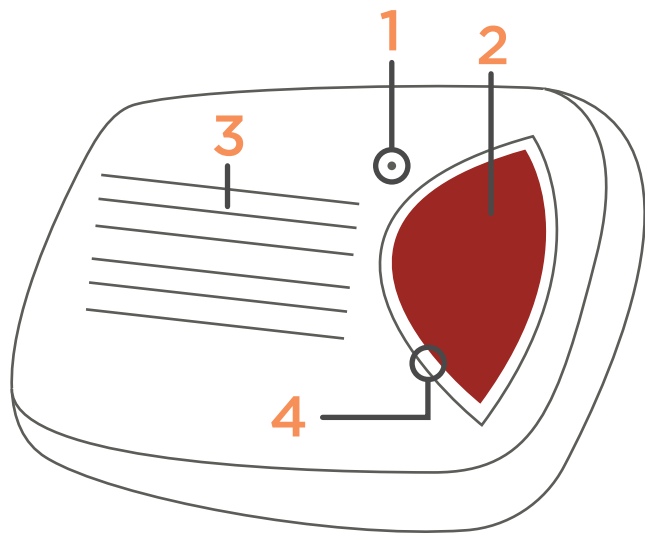

- 1 Microphone**
Used to talk to the Care Center
- 2 Help Button**
Press to call the Care Center (works just like the button)
- 3 Speaker**
Used to hear the Care Center
- 4 Status Light**
Displays the current operating condition
- 5 Cellular Antenna**
Provides integrated cellular service (no attachments necessary)
- 6 Programming Jack**
For service provider use only
- 7 Reset Button**
Used to Reset the unit
- 8 T/L Button**
The Test and Learn Button. Used in the programming and testing of the unit

Unit Facts & Features:

- ✦ The button can be worn around your neck or as a wristband
- ✦ The button has a range of up to 600 feet, (average 3-5 year battery life), and is water resistant
- ✦ Unit has its own internal cell modem. (Does not use clients cell plan)
- ✦ Care Center monitors unit's battery status
- ✦ 60+ hours of backup battery in case of AC power loss
- ✦ Optional fall detection pendant can be added
- ✦ Optional Medication monitoring can be added
- ✦ Status Light Colors:
 - **Solid Green-** Unit is properly connected
 - **Flashing Green-** AC power is not available, unit is running on back up battery
- **Steady Red-** Unit is trying to call the Care Center
- **Flashing Red-** unit does not have sufficient cell signal
- **Flashing Orange-** Unit does not have sufficient cell signal and is running on backup battery
- ✦ For additional questions, call tech support at **1-800-860-4230**