


Meet Your New MDS

MedReady Medication Dispenser System

MedReady was created to help manage your medication to always take the right dose at the right time.

Get to Know Your Unit

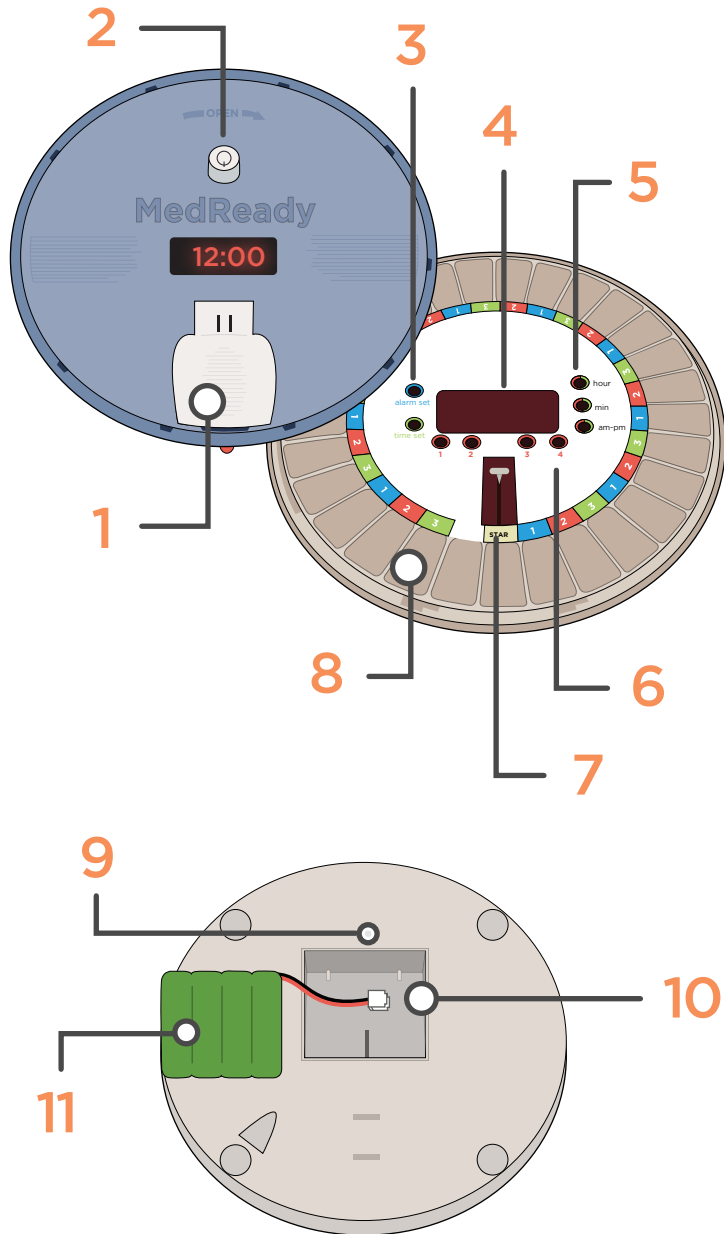

- 1 Medication Door**
Allows client to remove medication and shuts off audible alarm
- 2 Lock and Key**
Allows locking of the MedReady
- 3 Alarm Test & Time Set Buttons**
Alarm test changes the amount of time the audible alarm will go off
Time set will enable the local time to be set on the clock
- 4 Clock and Indicator Light**
Displays the local time and alarm times
- 5 Hour, Minute, AM/PM Buttons**
Changes the Hour, Minute and AM/PM for local time and alarm times. AM/PM will also test the device
- 6 Alarm Set Buttons (1,2,3,4)**
Set alarm times (1 - 4)
- 7 White T-Bar (Resets Alarm)**
Resets alarm - Bar must be pushed up toward time before placing lid on to work properly
- 8 Medication Tray**
Medication placeholder
- 9 Power Jack**
Connects unit to power outlet
- 10 Battery Jack**
Connects battery pack to unit
- 11 Battery Pack**
Back-up battery to unit

Unit Facts & Features:

- ✦ MedReady holds up to 28 doses
- ✦ Dispenses up to four times a day
- ✦ Locking lid prevents double doses
- ✦ For additional questions, call tech support at [1-800-860-4230](tel:1-800-860-4230)